





NextGENe[®]

Network Version Setup Guide

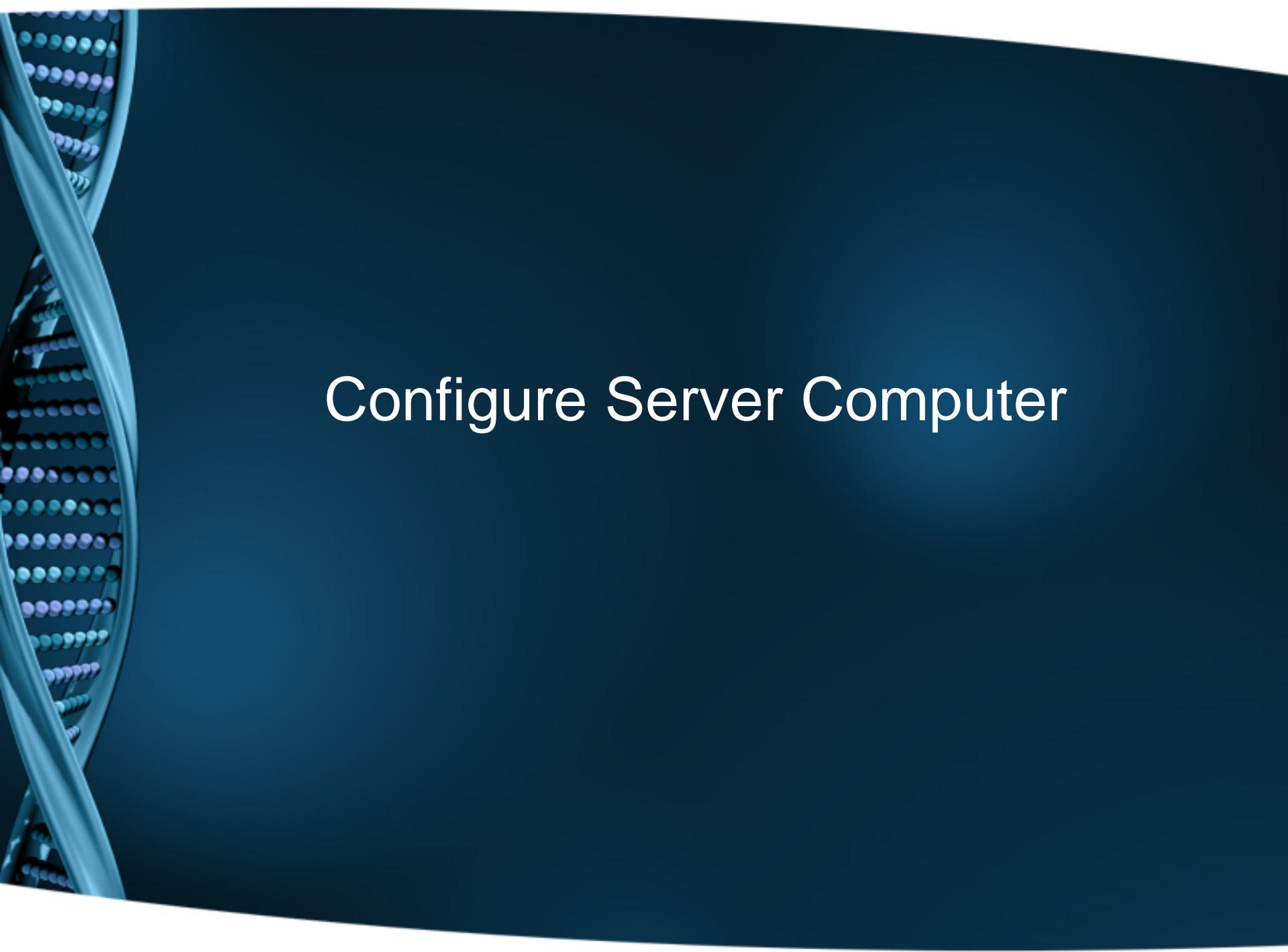


Configure Server Computer

1. Log on to the Windows Server as a Local Administrator. In order to meet User Account Control settings, you may also need to right-click on the NextGENe setup and choose Run As Administrator.
2. Choose to Install License Server Manager.
NOTE: This program monitors available licenses from clients that are running NextGENe. It should be installed on a Server or a computer that is always running. The License Server Manager can be installed on the same computer that will be running NextGENe.

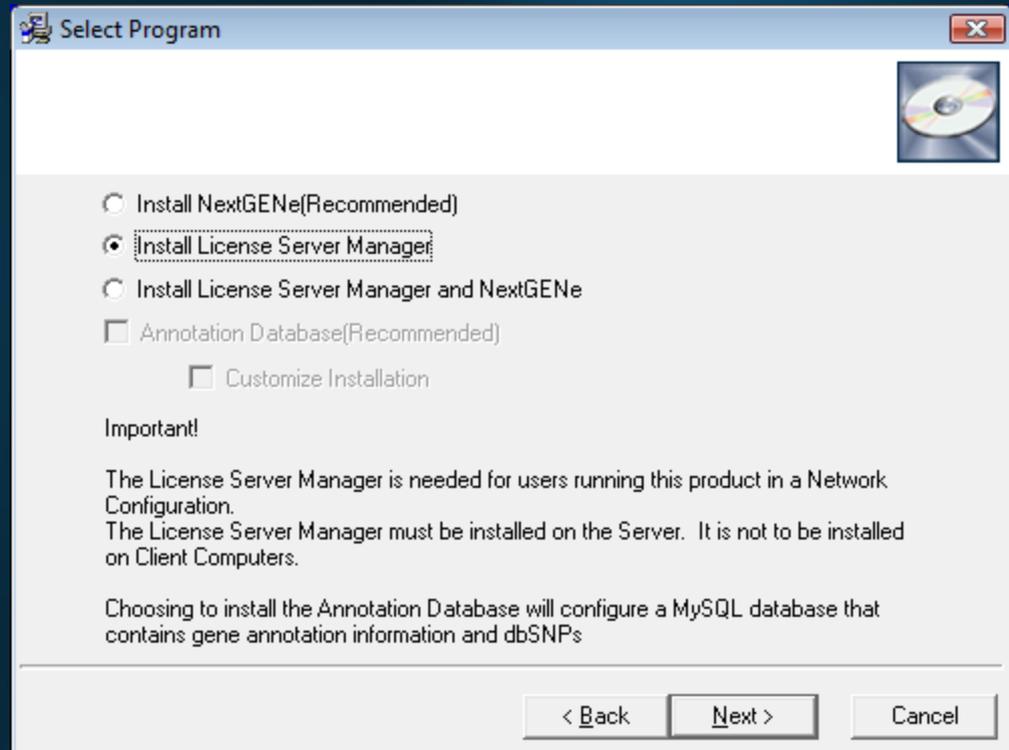


1. Click Next to Install License Server.
2. When complete, choose to Launch License Server.

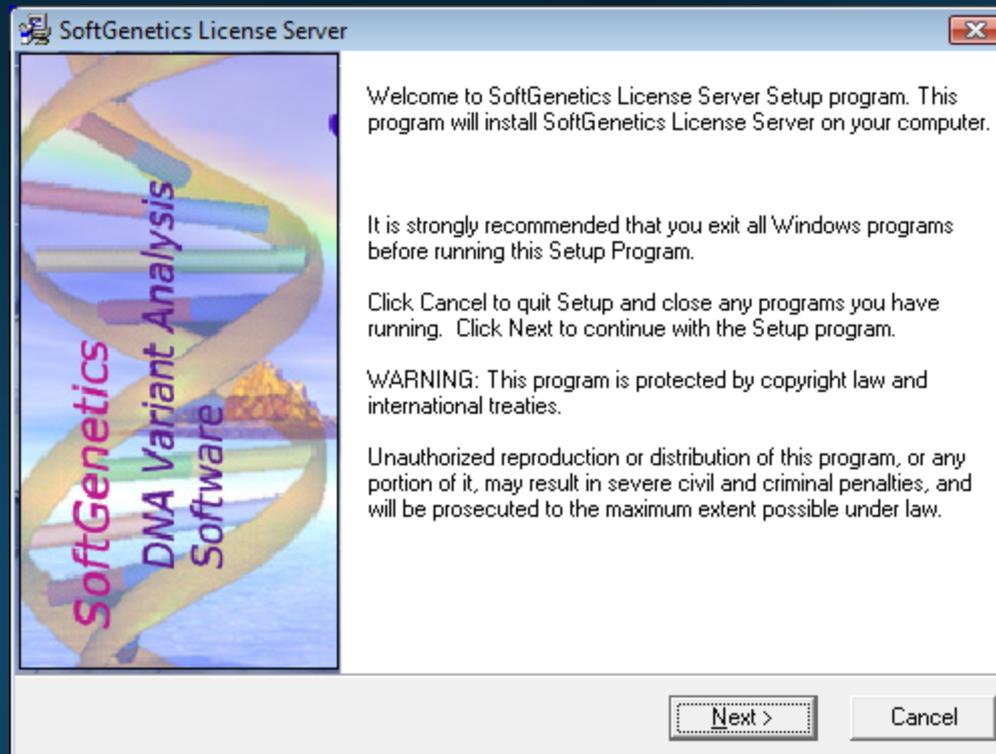



Register License Server Manager for NextGENe

1. Open License Server Manager  from the Icon Tray by clicking on LSM icon in System Tray
2. Register License Server for NextGENe by clicking on “Register” from the Help menu
3. Select “NextGENe” as Product
Input your Account, Password and Email and click Register
4. You will be prompted with a message to restart License Server Manager for the changes to take effect

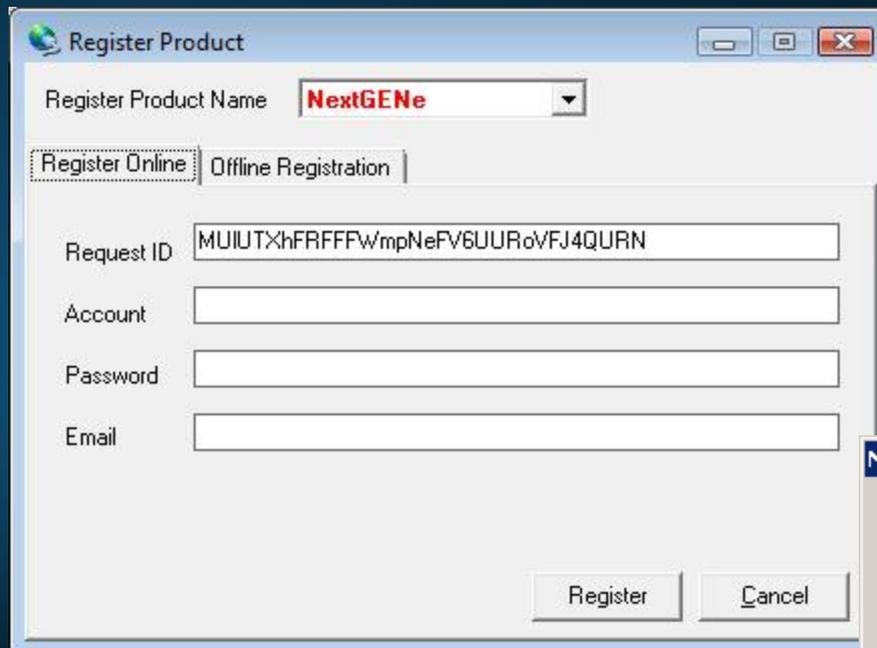

Register Product

Register Product Name: **NextGENe**


Register Online Offline Registration

Request ID: MUIUTXhFRFFFWmpNeFV6UURoVFJ4QURN

Account:


Password:

Email:

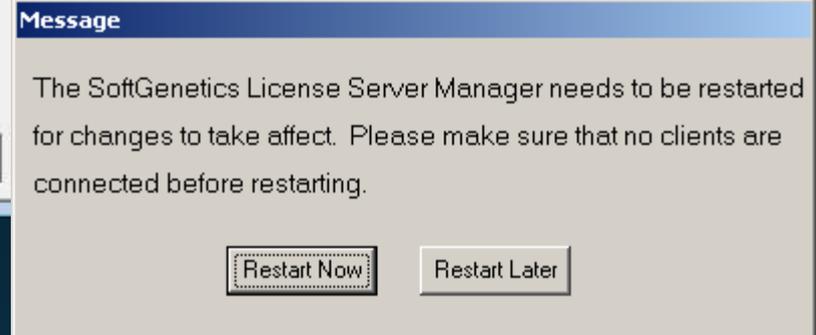


Message


The SoftGenetics License Server Manager needs to be restarted for changes to take affect. Please make sure that no clients are connected before restarting.

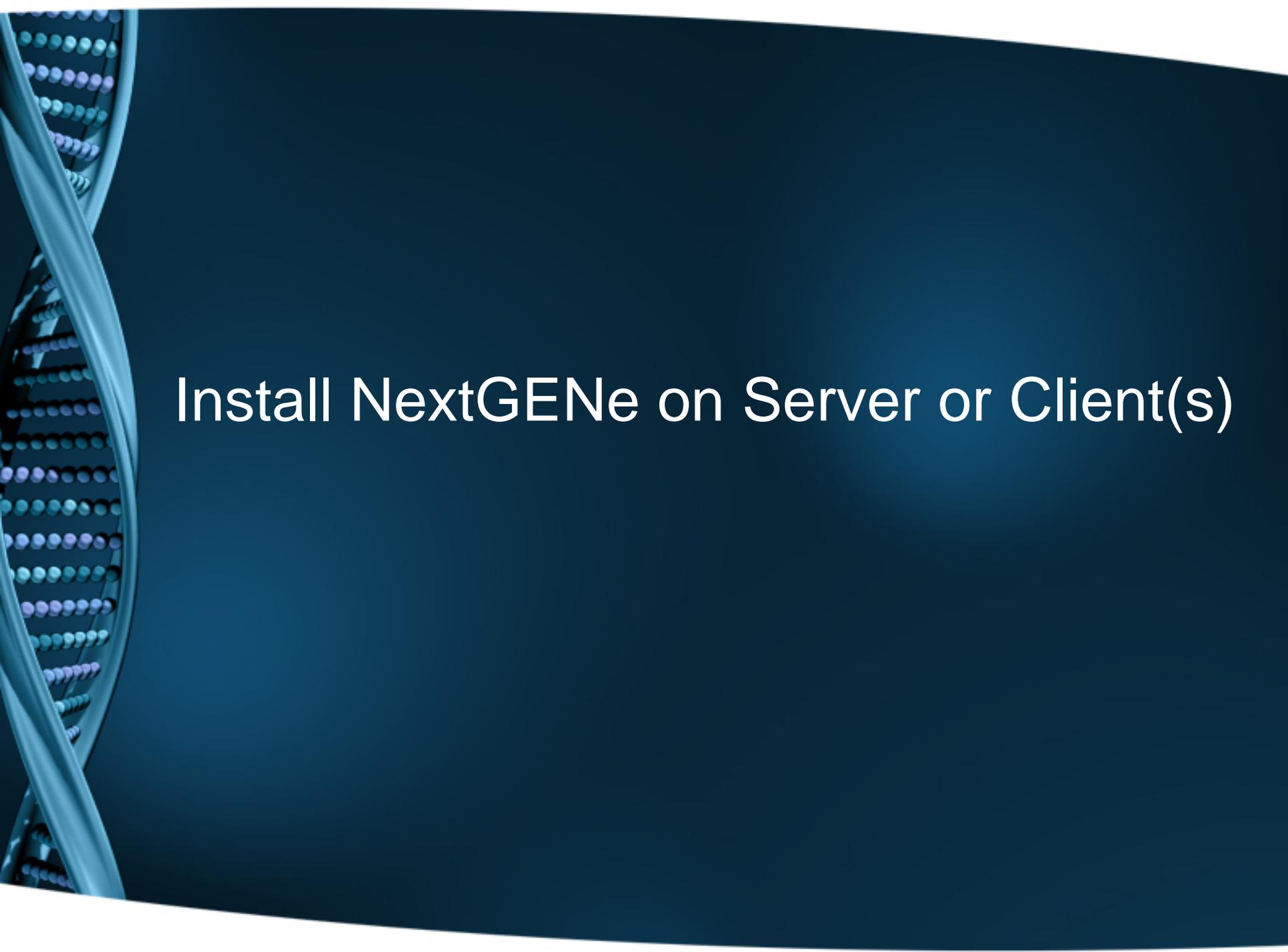

Install NextGENe on Server or Client(s)



1. Log on to the client computer as a Local Administrator. In order to meet User Account Control settings, you may also need to right-click on the NextGENe setup and choose Run As Administrator.
2. Choose to install NextGENe.

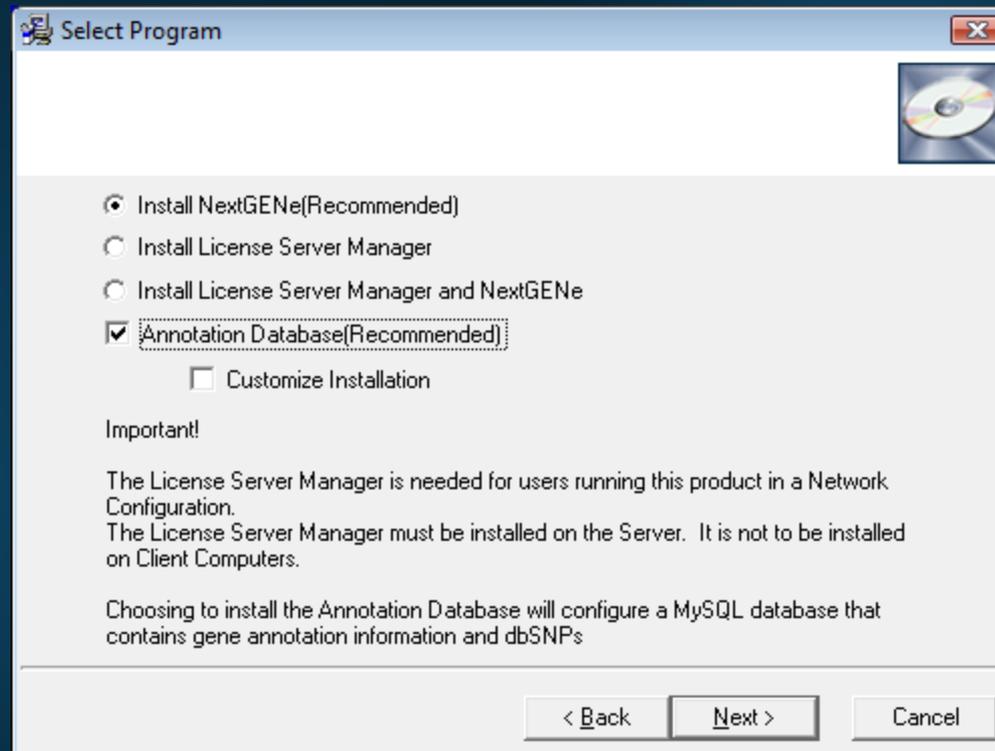

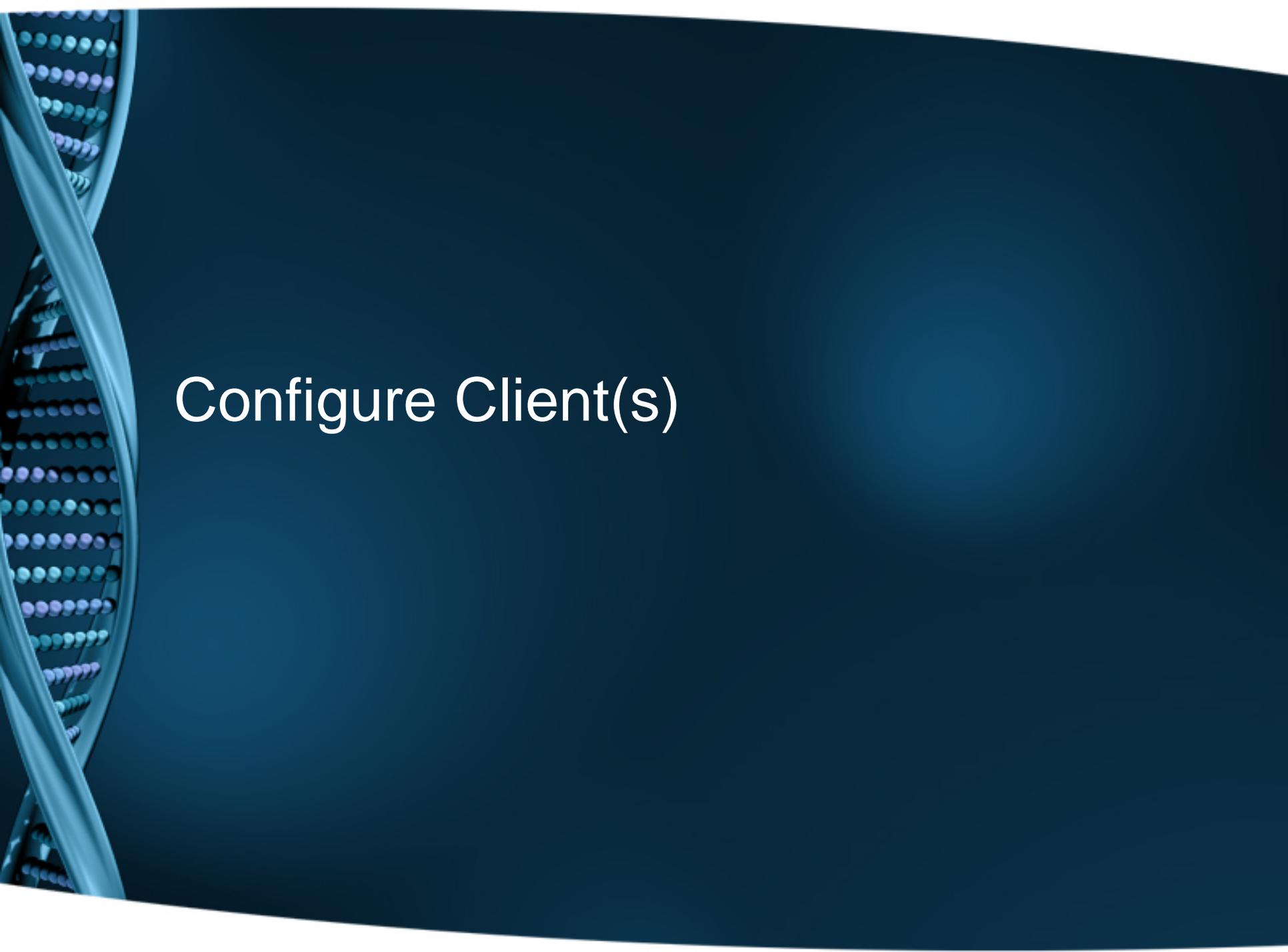

Configure Client(s)

1. Open NextGENe. “License is not configured” message will appear.
2. Click OK.
3. Choose Network type and click Configure.

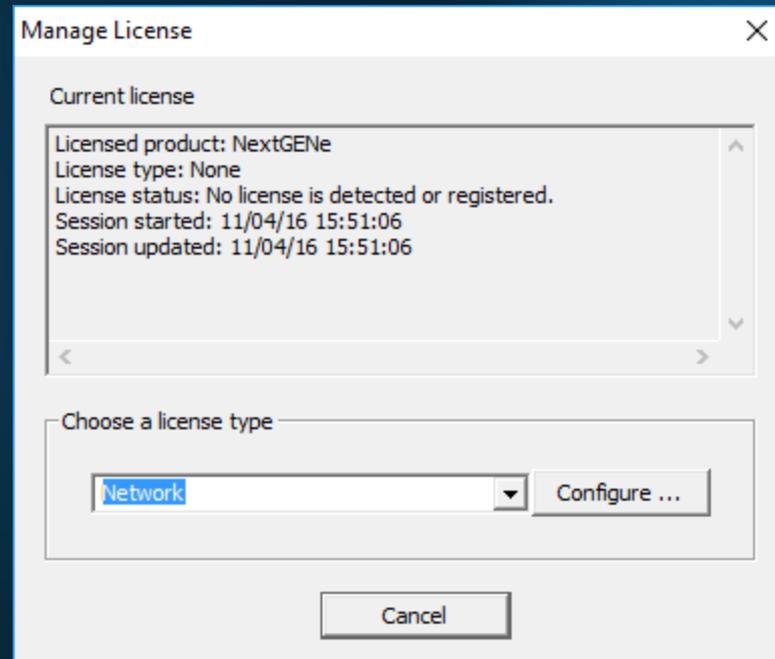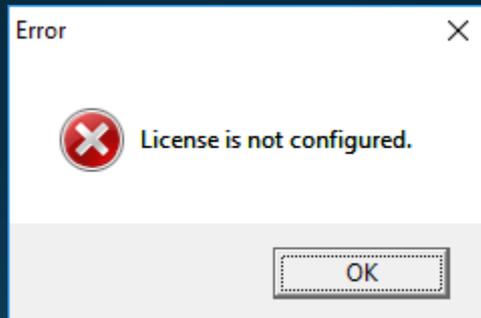

1. Input host name or IP Address.
2. Click Connect. Message will show “Succeeded”.
3. Click OK.
4. NextGENe will automatically open if connection is properly established and a license is available.

