

NextGENe®


Network Upgrade Guide for License Server Manager



Update License Server Manager on Server Computer




- Log on to the Windows Server as a Local Administrator. In order to meet User Account Control settings, you may also need to right-click on the NextGENe setup and choose Run As Administrator.
- 2. Choose to Install License Server Manager.

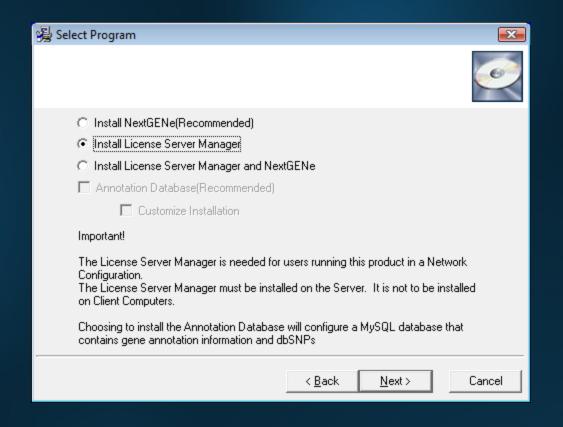



- Click Next to Install License Server.
- 2. When complete, choose to Launch License Server.

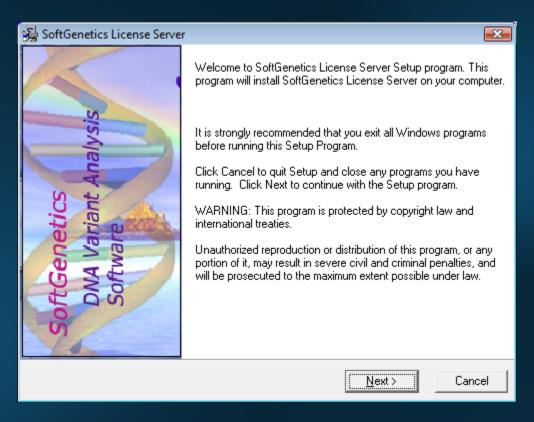



Register License Server Manager for NextGENe

- This is needed only when the valid period is expired.
- If the previous version of NextGENe was running on clients without any problems, then this step is not necessary



- 1. Open License Server Manages from the Icon Tray by clicking on LSM icon in System Tray
- 2. Register License Server for NextGENe by clicking on "Register" from the Help menu
- 3. Select "NextGENe" as Product Input your Account, Password and Email and click Register
- 4. You will be prompted with a message to restart License Server Manager for the changes to take effect


Register Product	
Register Product Name NextGENe ▼	
Request ID MUIUTXhFRFFFWmpNeFV6UURoVFJ4 Account Password Email	QURN
F	The SoftGenetics License Server Manager needs to be restarted for changes to take affect. Please make sure that no clients are connected before restarting. Restart Now Restart Later



Install NextGENe on Server or Client(s)



- Log on to the client computer as a Local Administrator. In order to meet User Account Control settings, you may also need to right-click on the NextGENe setup and choose Run As Administrator.
- Choose to install NextGENe.

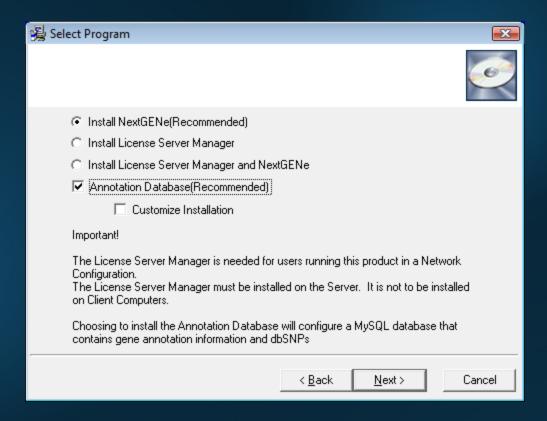




NextGENe will now work on previously configured client computers

 To configure new client computers, continue with next steps



- 1. Open NextGENe. "License is not configured" message will appear.
- 2. Click OK.
- 3. Choose Network type and click Configure.



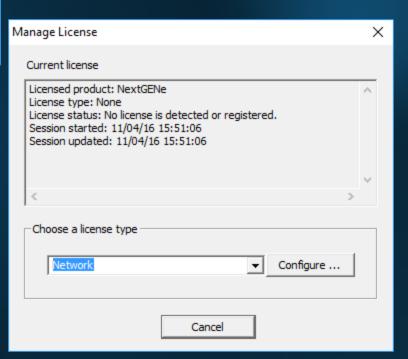



- Input host name or IP Address.
- 2. Click Connect. Message will show "Succeeded".
- 3. Click OK.
- 4. NextGENe will automatically open if connection is properly established and a license is available.

